



# PROTEOMIC

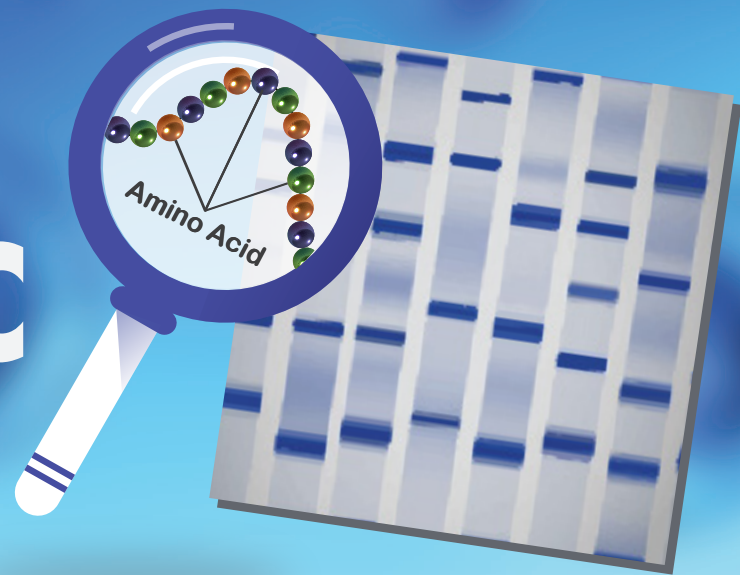

## PROTEIN IDENTIFICATION SERVICE

Proteomic protein identification service refers to a specialized service offered by Genetika Science to identify proteins present in a biological sample. Proteomics is the study of proteins and their functions within an organism or a biological system.

### WORKFLOW

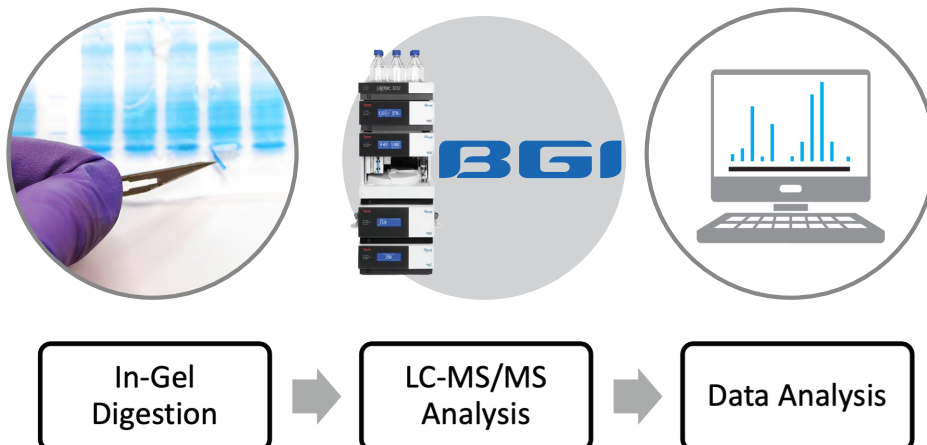

Gel band identification uses liquid chromatography with tandem mass spectrometry (LC-MS/MS) method to identify a protein band in sodium dodecyl sulfate polyacrylamide (SDS-PAGE).

This method combines the physical separation capability of liquid chromatography with the mass analysis capabilities of mass spectrometry. Furthermore, this method utilizes tandem mass spectrometry where two or more mass analyzers are coupled together increasing their abilities to analyze protein sample.

## Key Deliverables for Gel Band Protein Identification:

- **Summary report including all methods and data analysis.**
- **Result files and protein sequences in fasta format.**
- **Bioinformatic analysis:**
  - Identified peptide and identified protein list.
  - Functional Annotation (GO, KOG, KEGG).
- **Raw files available upon request.**

## SAMPLE REQUIREMENT

<b>Platform</b>	: Ultimate 3000 UHPLC & Q Exactive HF-X Quadrupole Orbitrap MS System.
<b>Sample type</b>	: Protein solution.
<b>Protein amount:</b>	: $\geq 400\mu\text{L}$ protein solution ( $\geq 0.5\ \mu\text{g}/\mu\text{L}$ ).
<b>Minimum sample amount</b>	: $\geq 200\mu\text{L}$ protein solution ( $\geq 0.5\ \mu\text{g}/\mu\text{L}$ ).
<b>Storage condition</b>	: $-80\text{C}$ or dry ice (for transport only).

### Note:

- Samples does not contain loading dye or SDS.
- Samples that require enrichment are not accepted.



**PT. Genetika Science Indonesia**

Green Lake City, Rukan Great Wall Blok C No. 19-21,  
Kel. Petir, Kec. Cipondoh, Kota Tangerang,  
Banten 15147, Indonesia.

SCAN HERE

